


IntelliStudio®




DermaGraphix®
body mapping solution

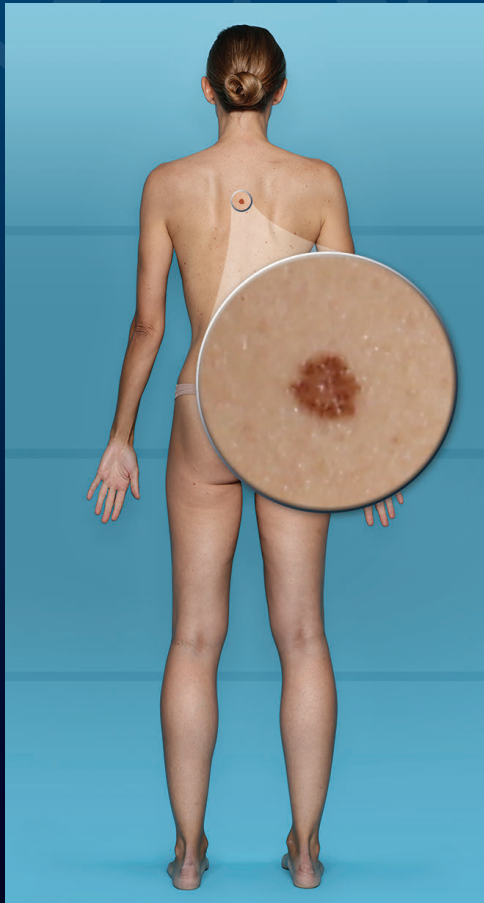

*50MP camera allows
for high magnification
of the smallest skin
features.*



*Close-ups with
non-polarized (left) and
cross-polarized lighting.*

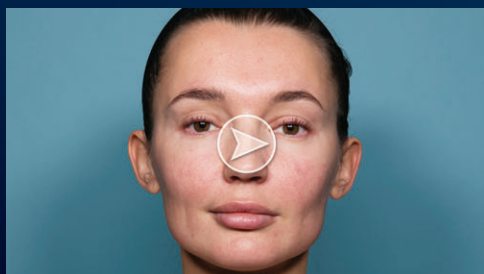

*Video with LED lighting
and record, playback
and zoom capabilities.*



IMAGING EXCELLENCE FROM
 **CANFIELD**




DermaGraphix

body mapping solution

Fast and repeatable total body photography with IntelliStudio and DermaGraphix body mapping software.

- Fully automated image capture for consistent clinical quality images.
- Templates control the camera settings and guide patient poses.
- Motorized programmable lift ensures repeatable camera positions and lighting for consistent, reproducible images.

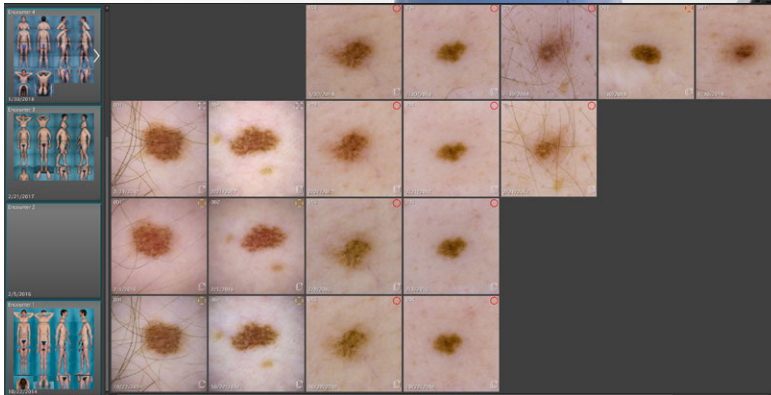


The DermaGraphix Dashboard organizes patient images for optimal review and assessment.

- View all captures of the same lesion in a single column; all lesions captured during the same visit in a single row.
- Secure HL7 compatible image management system.



Tag dermoscopic images to overview photos for efficient exams and accurate follow-ups.

- All images associated with a lesion are linked to a unique number and tagged to the body map.
- Selectable tagging routines provided to accommodate workflow preferences.
- All photos, notes and attributes associated with a lesion are available for instant review and comparison.


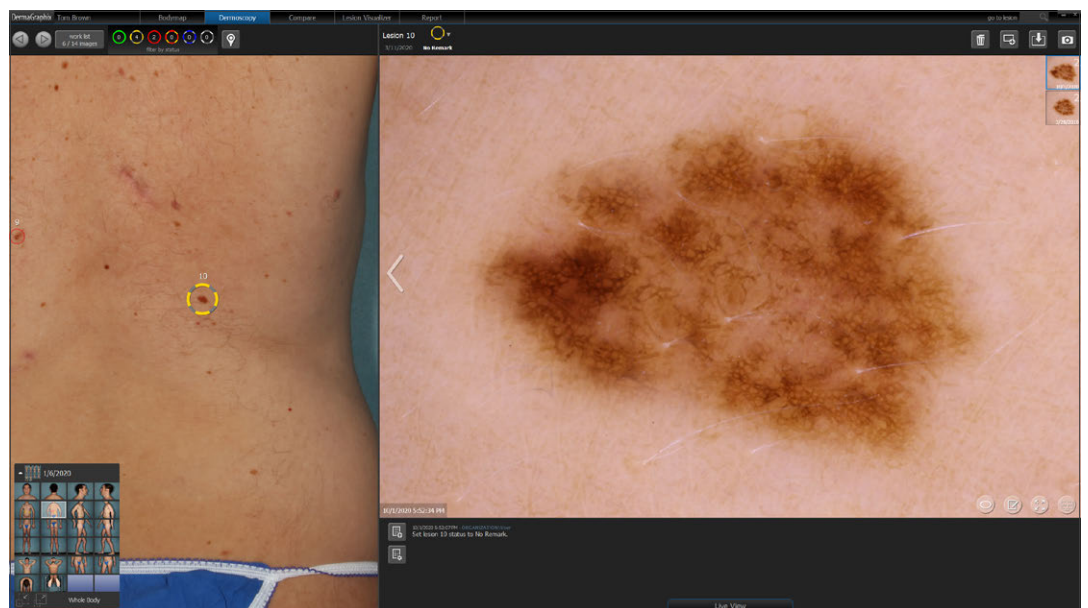


Smart Compare and Lesion Visualizer enhance change detection

Smart Compare uses Canfield's Markerless Tracking™ to identify corresponding regions between sets of images and bring them into ideal alignment, with each lesion matched across timepoints. This enables effortless visual comparison of two time points by effectively correcting for patient pose variation.

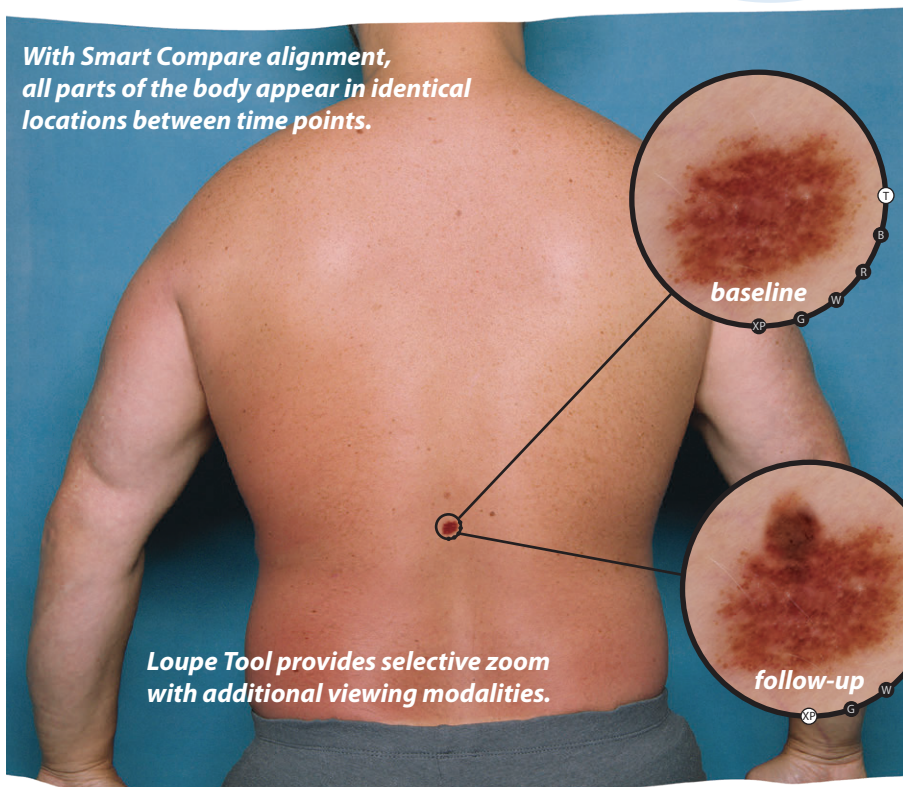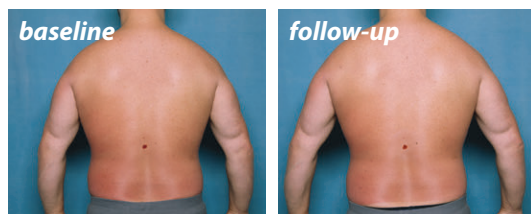

Lesion Visualizer* automatically detects, analyzes, and tracks lesions on bodymap overview images to reveal evolving lesions.

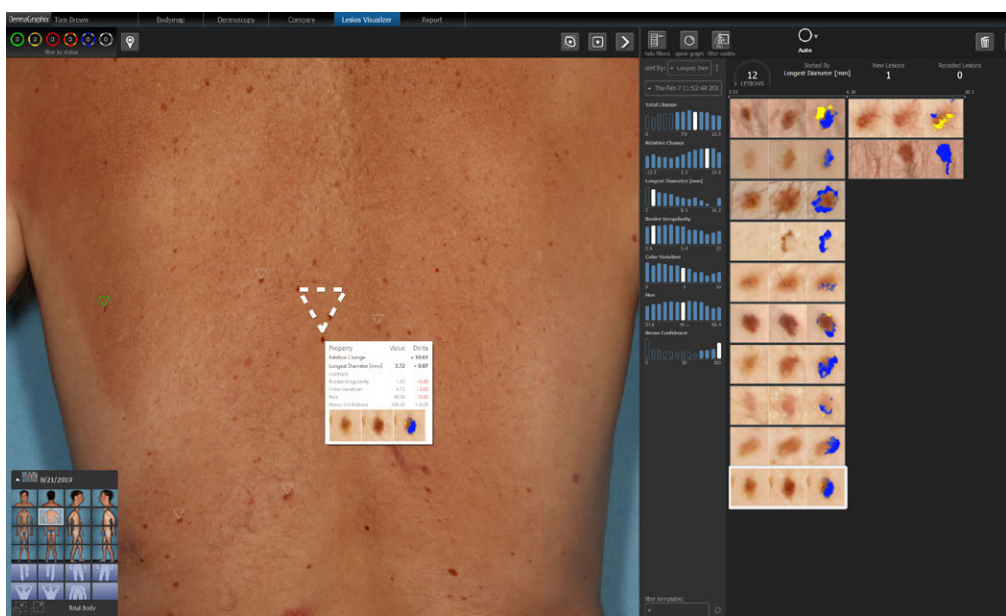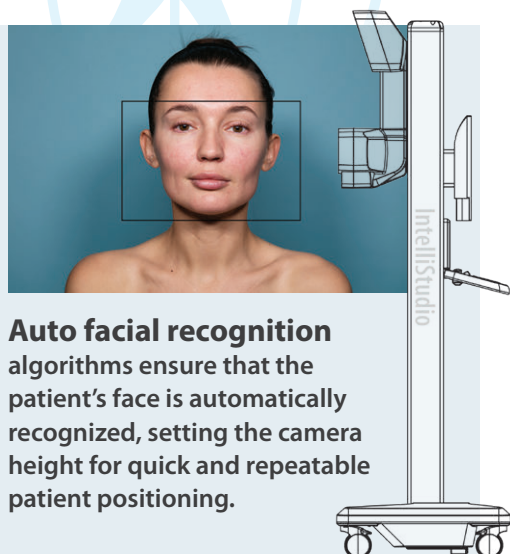

*Lesion Visualizer is for research purposes only.

IntelliStudio®

Studio-quality xenon lighting Using a 50MP camera and studio-quality lighting with smart automation, non-polarized and cross-polarized lighting is calibrated and balanced for skin imaging, revealing skin features while maintaining aesthetics and reducing shadows.


Motorized programmable lift With automated sequencing to guide the user through captures using standardized photo templates, the motorized programmable lift ensures repeatable camera positions and lighting for consistent, reproducible images. Integrated DermaGraphix® software enables quick capture of body map photos.

Video capabilities with LED lighting IntelliStudio provides video capabilities, including record, playback and zoom capabilities.

MatchPose® overlay ensures registered before-and-after images.

Fully automated IntelliStudio simplifies image capture with an automated software process that includes a selection of standard templates. With IntelliStudio, there are no special photography skills needed. Templates control the camera settings and guide patient poses, ensuring consistent clinical quality images every time.

Programmable zoom allows you to set and save specific program settings, enabling consistent zoom levels for accurate, repeatable close-up photos.

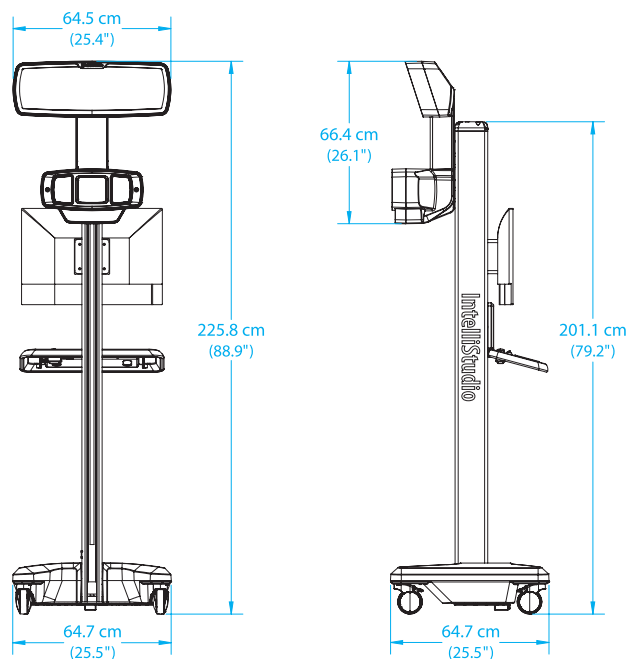

Powerful search capability with customizable attributes for patients, lesions and images. Search by doctor, procedure, diagnosis, date ranges and other attributes. Audit trail and stratified rights capable.

DermaGraphix Viewer empowers patients to check for changes between office visits. The simple, intuitive user interface includes both overview and close-up images. Data is stored on a secure, password protected and encrypted USB flash drive.

CANFIELD dermatoscopes

Canfield dermatoscopes are compatible with DermaGraphix software so you can easily link each close up image to its corresponding tagged lesion in the body map.

VISIONED™ D200evo

A high-resolution digital dermatoscope with quality, flexibility and efficiency all in a compact, hand-held system.

- 15–200x magnification
- Polarized and non-polarized lighting
- Tag and capture using the touch screen and keypad for computer-free dermoscopy photos
- Capture both clinical (overview) and dermoscopy photos



VEOS® DS3 Canfield's exclusive optics and lighting combined with wireless image capture.

- Sequentially capture cross-polarized and non-polarized images
- Capture both clinical (overview) and dermoscopy photos
- Contact and non-contact imaging
- 10x magnification
- Integrated iPod® Touch for touch screen operation



VEOS SLR captures consistent, reproducible dermoscopic images.

- 22mm field of view and high resolution capture.
- Sequentially capture cross-polarized and non-polarized images in either contact or non-contact modes.



IMAGING EXCELLENCE FROM



www.canfieldsci.com / info@canfieldsci.com / phone +1.973.434.1201 / (USA) 800.815.4330

3D SOLUTIONS • FACIAL IMAGING & ANALYSIS • IMAGING SOFTWARE • PHOTOGRAPHY • RESEARCH SYSTEMS & SERVICES • TRAINING

IntelliStudio, DermaGraphix, MatchPose, Markerless Tracking, VEOS, and VISIONED are trademarks of Canfield Scientific, Inc.
iPod is a trademark of Apple Inc., registered in the U.S. and other countries.

2009-25