



EVIS EXERA III VIDEOCOLONOSCOPE

PCF-H190TL/I

High-definition colonoscope facilitating endoscopic treatment through excellent maneuverability



PCF-H190TL/I



Main Features

Short bending for better access to lesions





With its HDTV image quality and

brighter NBI, the PCF-H190TL/I

contributes to higher treatment

diagnosis facilitates appropriate

delineation of lesions.

accuracy. High confidence optical


treatment decisions and accurate

PCF-H190TL/I

HDTV and **NBI**

White Light

Appropriate channel layout

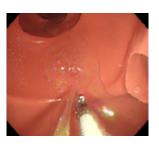

The layout allows for smooth operation of both the water-jet function and any treatment such as stopping bleeding.

Improved insertability



PCF-H190TL/I

The PCF-H190TL/I incorporates unique proprietary Olympus technologies, realizing easier scope insertion with minimized patient discomfort.

Specifications Field of view (Maximum diagonal angle) 140° Minimum visible distance 4 mm from the distal end Direction of view 0° (Forward viewing) Instrument Direction from which endotherapy **Optical System** Channel accessories enter and exit the Depth of field 2-100 mm endoscopic image Direction of water jet in field Distal end outer diameter 9.8 mm Water Jet of view Up 210º/Down 180º/ Distal end enlarged Bending Section Angulation range Right 160°/Left 160° Up Objective Lens Air/Water Nozzle **Total Length** L: 2005 mm, l: 1655 mm Insertion Section Right Left Compatible Light-guide Lens Video System Center Olympus CV-190 Instrument Channel Outlet EVIS EXERA III Auxiliary Water Channel Xenon Light Source Olympus CLV-190 System Down Insertion tube outer diameter 10.5 mm Working length L: 1680 mm, l: 1330 mm Instrument Channel Channel inner diameter 3.2 mm

Manufactured by Aomori Olympus Co., Ltd, 2-248-1, Okkonoki Kuroishi, Aomori, 036-0357 Japan. Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer. Olympus is a registered trademark of Olympus Corporation, Olympus America Inc., and/or their affiliates. I Medical devices listed may not be available for sale in all countries.



OLYMPUS AMERICA INC. 3500 Corporate Parkway, PO BOX 610, Center Valley, PA 18034 For more information, contact your Olympus sales representative, or call 800-848-9024. www.medical.olympusamerica.com

©2020 Olympus America Inc. All rights reserved. Printed in the USA OAIGI0820BRO36282